
Gebrauchsanweisung | Instructions for use |
Mode d'emploi | Modo de empleo

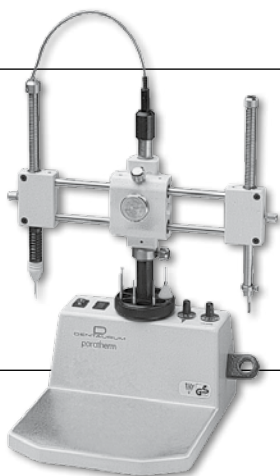Paratherm

REF 094-270-00

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Druckdatum: 28.10.15

Dear customer

Thank you for choosing a quality product from Dentaurum.

It is essential to read these instructions carefully and adhere to them to ensure safe, efficient use and ensure that you and your patients gain full benefit. Instructions for use cannot describe every eventuality and possible application. In case of questions or ideas, please contact your local representative.

As our products are regularly upgraded, we recommend that you always carefully read the current instructions for use supplied with the product and stored in the internet at www.dentaurum.com, even though you frequently use the same product.

1. Manufacturer

Dentaurum GmbH & Co. KG | Turnstraße 31 | 75228 Ispringen | Germany

2. With reference to our quality standards

Dentaurum ensures a faultless quality of the products manufactured by us. These recommendations are based upon our own experiences. The user himself is responsible for the processing of the products. Responsibility for failures cannot be taken, as we have no influence on the processing on site.

3. Safety instructions

Caution

Read through the instructions for installing and operating the unit carefully before switching it on.

3.1 Disposal

When a Dentaurum product is finally no longer in use the disposal regulations of that particular country apply. Dentaurum or the dental trade is available to answer questions regarding the correct disposal of any specific product.

Note:

This symbol indicates, products marked with it should not be disposed off together with garbage of households. The legislator disallows commercial customers to return electronic waste via municipal collection points. Further information is available through Dentaurum or the dental trade.

3.2 Delivery

	Description		piece
	Paratherm unit with surveying / blocking-out arm		1
Ⓔ	Thermal blocking-out instruments 0°, 2°		1
Ⓓ	Undercut gauges 0.25 / 0.35 / 0.50 / 0.75 mm	0,25 mm	4
Ⓐ	Analyzing rod		1
Ⓑ	Graphite holder		1
Ⓒ	Replacement graphite markers		10
	Plastic cover		1
	Fuse T 0.25 A		1
	Instructions for use with EEC-Declaration of conformity		10

4. Areas of use

Cast partial denture, combination techniques, conical crowns, positioning of snap attachments and attachments.


5. Installation

Set the necessary voltage on the underside of the unit (110 / 220 / 240 V). Then connect the unit to the mains supply. Stand the Paratherm on a level working area.

6. Operating the Paratherm

Mount the master model on the model table already available or supplied with the unit. Loosen the knurled screw ② and raise the sliding rods to the necessary level. Fasten with retainer ring ①. Grip the instrument eg. analyzing rod in the surveying handpiece ④ and push the selected blocking-out instrument into the thermal blocking-out handpiece ⑤.

Place the model table ⑥ in the most convenient position on the work surface ⑦. Loosen the locking screws ③, so that the arms may be moved freely. Switch on the thermal blocking-out handpiece ⑧ and set the desired temperature (maximum 160 °C / 320 °F) on the controller ⑨.

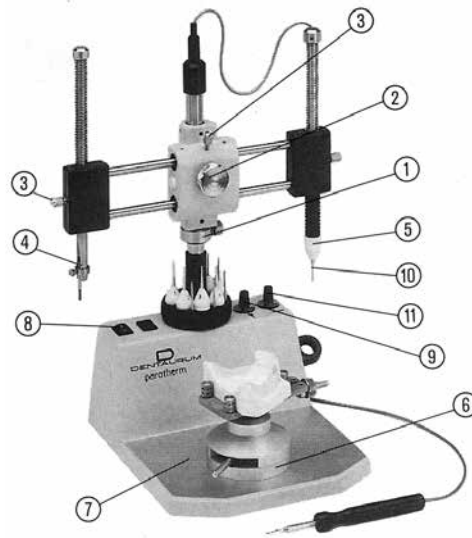Caution: Only change the blocking-out instruments ⑩ whilst they are still warm.

6.1 Electric wax knife

The knife is connected by means of the socket on the rear side of the unit. Only operate the wax-knife with an instrument in the place.

Set the desired temperature on the potentiometer ⑪. Maximum temperature 280 °C / 536 °F.

Caution! Only change the knife whilst it is still warm.

6.2 Locking the surveying arm (Positioning attachments)

Rotate the surveying handpiece ④ forwards of the centre of the work surface and lock into the position with screws ② ③.

Grip the wire aid in the surveying arm and fix the level using ③. With the model already clamped into position, guide the model table onto the wire aid. Correct the level.

When ordering spare parts please quote the unit number (type plate or rear side of the unit).

7. Spare parts

	Description		piece	REF	
Ⓔ	Thermal blocking-out instrument	0°	1	094-275-00	
		2°		094-277-00	
Ⓓ	Undercut gauge	0,25 mm	1	094-283-00	
		0,35 mm		094-284-00	
		0,50 mm		094-285-00	
		0,75 mm		094-286-00	
Ⓐ	Analyzing rod		1	094-287-00	
Ⓑ	Graphite holder		1	094-288-00	
Ⓒ	Replacement graphite markers		10	094-289-00	

8. Recommended accessories

(Not included in original kit REF 094-270-00)

	Description		piece	REF	
Ⓔ	Thermal blocking-out instrument	1°	1	094-276-00	
		3°		094-278-00	
		4°		094-279-00	
		6°		094-280-00	
		8°		094-281-00	

	Model table			086-340-00	
	Electric wax knife Handpiece, cable and 1 wax knife, holder for handpiece In order to regulate the temperature a transformer and potentiometer have been built into the Paratherm (REF 094-270-00)			094-272-00	
	Modelling knife, small			089-602-00	
	Thomas modelling point, fine			089-603-00	
	Thomas modelling point, medium			089-604-00	
	Orthodontic wax knife			089-605-00	
	Modelling spatula			089-606-00	

9.

EC-Declaration of Conformity

DENTAURUM GmbH & Co. KG
Turnstr. 31
75228 Ispringen

hereby declares that the design and construction of the laboratory equipment described below, including the version marketed by us, comply with the basic regulations governing safety and health as stated in the EC Guidelines. This declaration will become invalid if the laboratory equipment is modified or altered in any way without our prior consent.

Description of unit:	Paratherm	
REF:	094-270-00	
Start with unit No.:	32-080	
EC guidelines:	89/336/EEC 92/31/EEC	Electromagnetic compatibility Modified Guideline
	73/23/EEC	Electrical equipment used within certain voltage limits
	93/68/EEC guidelines	Modifications of the above
Applied unified standards:	EN 55011 EN 50082 EN 61010	

Date and manufacturers signature: 01.01.2010
Signatory:

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